


Tune Report

Operator Name ICPMS QQQ 8900
Acq/Data Batch D:\Agilent\ICPMH\1\DATA\2024\151024(1).b
Acq. Date-Time 10/15/2024 1:26:41 PM
Report Comment ---
Instrument Name G3665A SG23363340

[He]

Sensitivity

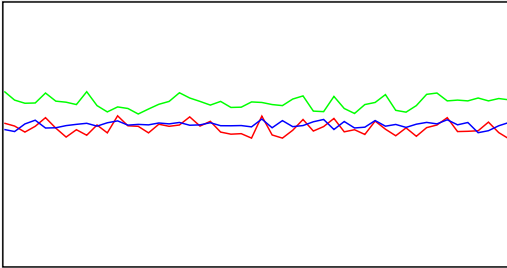

Mass	Range	Count	RSD%	Background
59	1000	524	4.559	
89	1000	621	3.625	
205	5000	2685	2.205	

Sampling Period [sec] 0.31
Integration Time [sec] 0.1

Oxide/Doubly Charged Ratio

Oxide 156 / 140 0.158 %
Doubly Charged 70 / 140 0.696 %

Tune Parameters

Plasma Parameters

Plasma Mode	---	Nebulizer Gas	0.75 L/min	Dilution Gas	0.40 L/min
RF Power	1600 W	Option Gas	0.0 %	Auxiliary Gas	0.90 L/min
RF Matching	1.60 V	Nebulizer Pump	0.10 rps	Plasma Gas	15.0 L/min
Sample Depth	10.0 mm	S/C Temp	2 °C		

Lens Parameters

Extract 1	0.0 V	Q1 Entrance	-50.0 V	Cell Exit	-60 V
Extract 2	-220.0 V	Q1 Exit	1.0 V	Deflect	-5.4 V
Omega Bias	-115 V	Cell Focus	-1.0 V	Plate Bias	-60 V
Omega Lens	5.6 V	Cell Entrance	-50 V		

Q1 Parameters

Q1 Mass Gain	127	Q1 Bias	-5.0 V
Q1 Mass Offset	125	Q1 Prefilter Bias	-6.5 V
Q1 Axis Gain	1.0002	Q1 Postfilter Bias	-9.0 V
Q1 Axis Offset	-0.07		

Q1 Ion Guide Parameters

SLS Factor	0.60	SLG Factor	0.90
------------	------	------------	------

Cell Parameters

Use Gas	Yes	3rd Gas Flow	0 %	Axial Acceleration	1.0 V
He Flow	5.0 mL/min	4th Gas Flow	0 %	OctP RF	180 V
H2 Flow	---	OctP Bias	-18.0 V	Energy Discrimination	3.0 V

Q2 Parameters

Q2 Mass Gain	125	Q2 Axis Gain	0.9998	Q2 Bias	-15.0 V
--------------	-----	--------------	--------	---------	---------

Tune Report

Q2 Mass Offset 126

Q2 Axis Offset 0.02

Q2 Ion Guide Parameters

SLS Factor 0.50

SLG Factor 1.00

Wait Time Offset

Wait Time Offset 0 msec

Hardware Settings

Torch

Torch H 0.3 mm

Torch V 0.1 mm

EM

Discriminator 2.7 mV

Analog HV 2183 V

Pulse HV 1153 V